

Charging options for ReSound hearing aids

ReSound offers rechargeable solutions with inductive charging and all-day battery life to keep your clients on the go.



Premium charger



Standard charger



Desktop charger



Custom charger

	Premium charger	Standard charger	Desktop charger	Custom charger
microRIE 60S ¹ 	8	8	8	-
miniRIE 60 ReSound Savi 	9	9	9	-
BTE 77 ² 	-	-	70	-
PBTE 88 ³ 	-	-	80	-
SPBTE 98 ReSound Enzo IA 	-	-	90	-
ITE/ITC ⁴ 	-	-	-	C

¹ **microRIE 60S:** ReSound Sensia, ReSound Vivia & ReSound Nexia

² **BTE 77:** ReSound Nexia

³ **PBTE 88:** ReSound Nexia

⁴ **ITE/ITC:** ReSound Sensia & ReSound Nexia

How to match ReSound hearing aids with the correct charger

Each charger type is offered with different inserts – known as the charging form factor (CFF) – which relate to a specific ReSound hearing aid collection or style.

The CFF is indicated with a number and is shown inside the charger, next to the serial number underneath the charger and on the charger package.





Charger features

	ReSound Sensia™, ReSound Vivia™ & ReSound Nexia™				ReSound Savi™		ReSound Enzo™ IA
	Premium charger	Standard charger	Desktop charger	Custom charger	Premium charger	Desktop charger	Desktop charger
	microRIE	microRIE	microRIE/ BTE/Power BTE	ITE/ITC	miniRIE	miniRIE	Super Power BTE
Inductive charging	•	•	•	•	•	•	•
Portable with onboard power bank	•				•		
Portable without power bank		•		•			
USB-C connector and cable	•	•	•	•	•	•	•
LED on the back	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply
LED on the front	Shows the hearing aids charging			Shows the hearing aids charging	Shows the hearing aids charging		

How to tell if the hearing aids are charging

Place the hearing aids in the charger and let the lights indicator on the hearing aids guide you. It is the same process whichever charger is used.

You know that charging has started when the hearing aid lights flash slowly.

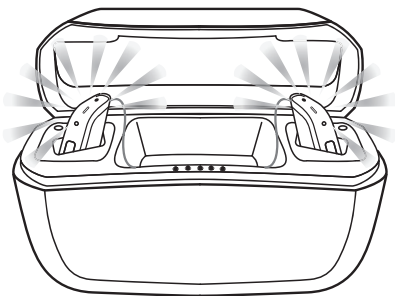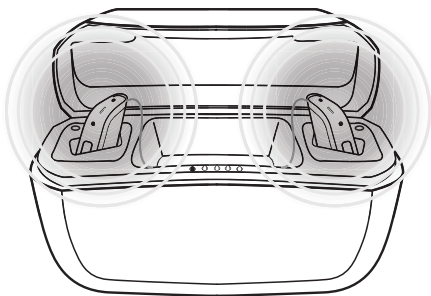

When fully charged, the hearing aid lights remain constant, not flashing.

The lights turn off when the hearing aids are removed from the charger.

Premium charger illustrated.
For more information, please see the individual user guide for each charger.

Worldwide Headquarters
GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: +45 4575 1111
gnhearing.com
CVR no. 55082715

United Kingdom
GN Hearing UK Ltd.
Unit 13
Talisman Business Centre
Bicester
Oxon OX26 6HR
Tel: +44 1869 352800
gnhearing.com

Australia
GN Hearing Australia Pty Ltd
Suite 1001, Level 1
1 Julius Avenue
North Ryde, NSW 2113
Australia
Tel.: (free) 1800 658 955
www.gnhearing.com/en-au/

New Zealand
GN Hearing New Zealand Limited
4 Fred Thomas Drive
Takapuna, Auckland, 0622
Tel.: (free) 0800 900 126
www.gnhearing.com/en-nz/